

# Aesculap Neurosurgery

## Yasargil Aneurysm Clip System

T-Bar Clips



Another essential tool for your  
Yasargil Aneurysm Clip tool box


# T-Bar Clip

## Aneurysm Clip with a unique design

The blades of the T-Bar Clip are oriented in a 90° (respectively 45°) angle to the axis of the spring plane resulting in a T-like clip form. The unique shape opens up new options for the treatment of complicated cerebral aneurysms. The spring type is identical to the existing and proven standard Yasargil Aneurysm Clip with a 5 mm diameter fenestration.

## Indications

- For clipping of cerebral aneurysms.
- T-Bar Clips can replace the use of multiple stacked fenestrated clips.
- For reconstruction of difficult to treat aneurysms, for example
  - complicated internal carotid artery aneurysms
  - fusiform aneurysms.

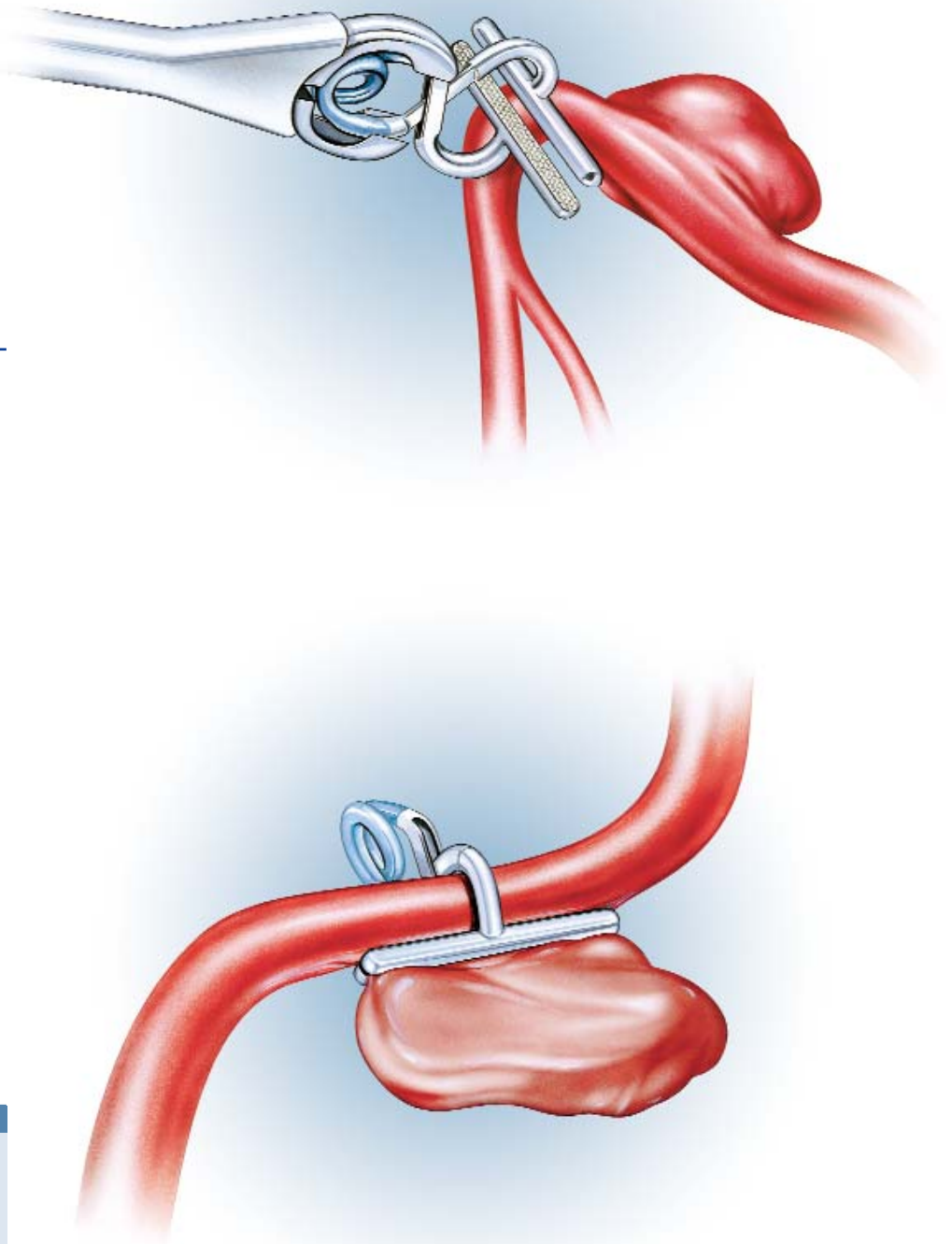

Developed in co-operation with:  
Dr. Thomas A. Kopitnik, M.D.  
Central Wyoming Neurosurgery  
Casper, Wyoming, USA

## Product Range

For **Permanent** implantation.

**Standard** size spring.

**Fenestrated** - fenestration diameter is 5 mm.



### T-Bar Clip Titanium

Titanium alloy acc. ISO 5832-3

Art. No.	Blade Length	Angle	Max. Opening	Closing Force		
FT852T	5 mm	45°	5.5 mm	150 g / 1.5 N		
FT853T	9 mm	45°	5.5 mm	150 g / 1.5 N		
FT854T	5 mm	90°	5.5 mm	150 g / 1.5 N		
FT855T	9 mm	90°	5.5 mm	150 g / 1.5 N		
FT856T	13 mm	90°	5.5 mm	180 g / 1.8 N		

### T-Bar Clip Phynox

Cobalt-base alloy acc. ISO 5832-7

Art. No.	Blade Length	Angle	Max. Opening	Closing Force		
FE852K	5 mm	45°	5.5 mm	150 g / 1.5 N		
FE853K	9 mm	45°	5.5 mm	150 g / 1.5 N		
FE854K	5 mm	90°	5.5 mm	150 g / 1.5 N		
FE855K	9 mm	90°	5.5 mm	150 g / 1.5 N		
FE856K	13 mm	90°	5.5 mm	180 g / 1.8 N		



AESCULAP®

**B | BRAUN**  
SHARING EXPERTISE

**Aesculap AG & Co. KG**

Am Aesculap-Platz  
78532 Tuttlingen  
Germany

Phone +49 7461 95-0  
Fax +49 7461 95-2600

[www.aesculap.de](http://www.aesculap.de)

All rights reserved. Technical alterations are possible. This leaflet may be used for no other purposes than offering, buying and selling of our products. No part may be copied or reproduced in any form. In the case of misuse we retain the rights to recall our catalogues and pricelists and to take legal actions.